





ARX CLAVICLE SYSTEM

2.5 / 3.5 Clavicle Plate

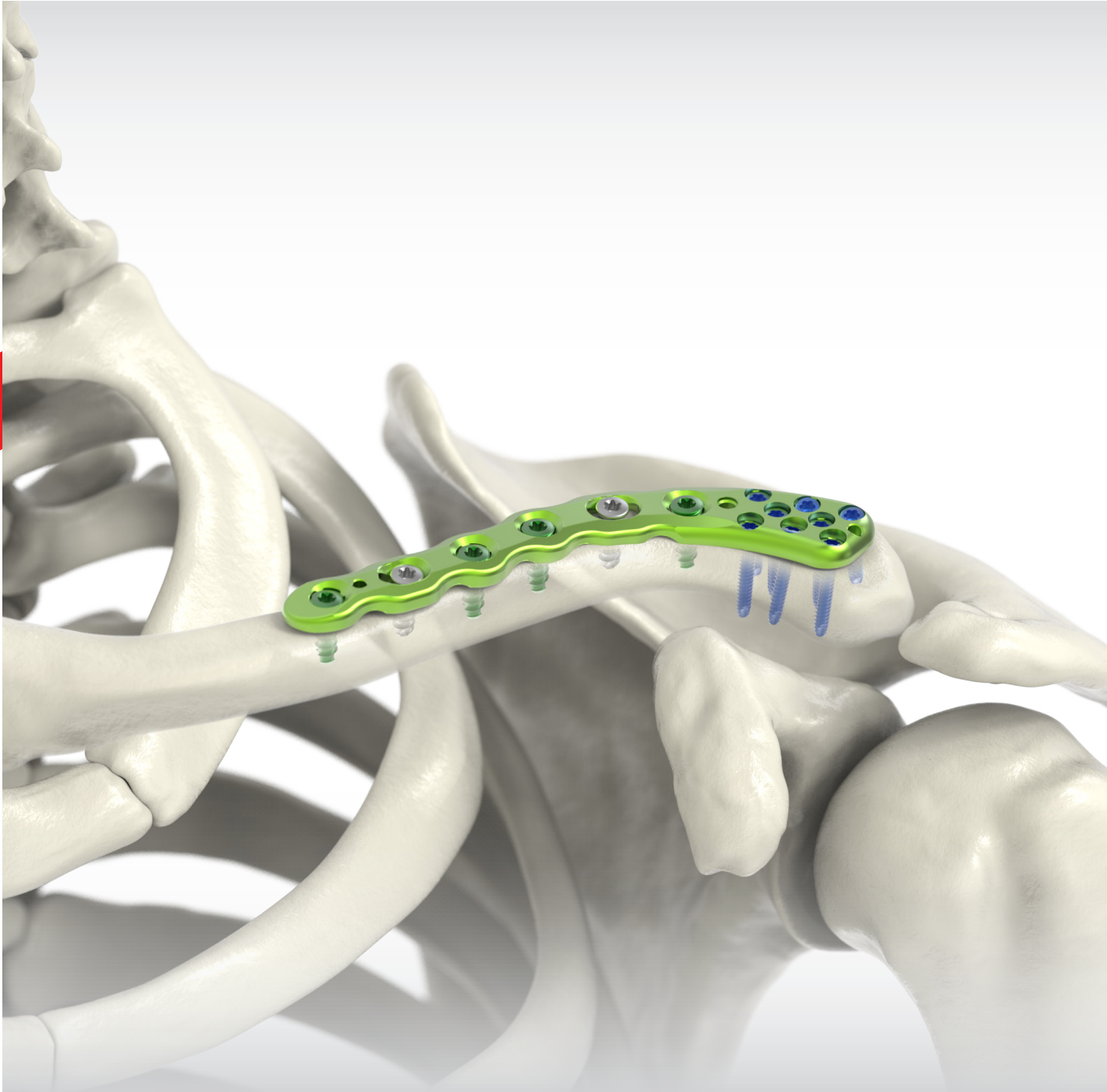

Contents

3	Indication
4	Overview
8	Ordering Information
8	- Clavicle Superior Lateral Plates
10	- Clavicle Superior Midshaft Plates
12	- Clavicle Hook Plates
14	- ARIX Clavicle System Screws
15	- ARIX Clavicle System Instruments
18	Set Configuration Suggestions

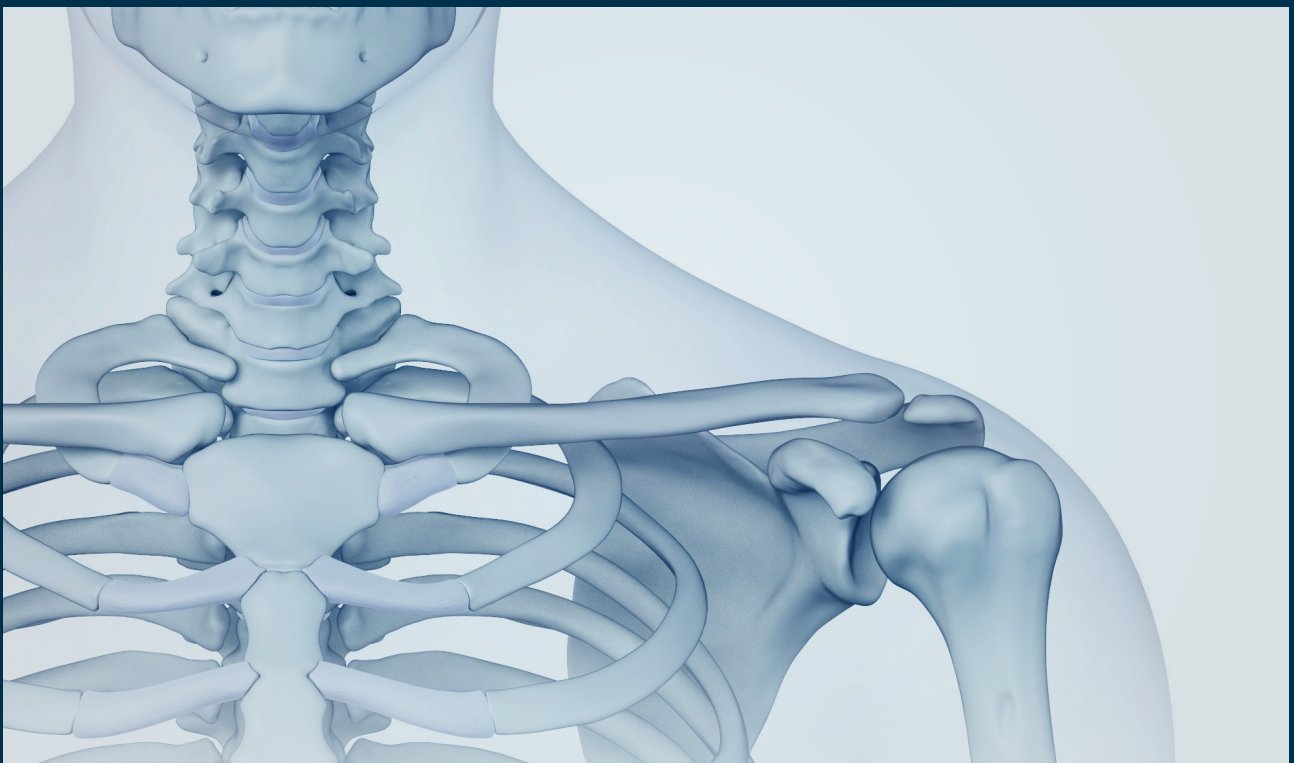


The ARIX Clavicle System is indicated for the fixation of single, segmental and comminuted fractures, osteotomies, mal-unions, and non-unions of the clavicle.

ARIX Clavicle Plate

The ARIX Clavicle System is a locking plate system for clavicle fractures. The Clavicle Plate provides an optimized fit for the clavicle bone. The ARIX Clavicle System also offers a variety of screws to help provide flexible surgical technique for users.

General Features

- Anatomically Pre-contoured
- Variable Locking up to $\pm 15^\circ$
- Left & Right Plate Variations for anatomical fitting
- Plate Thickness 3.5 mm
- $\varnothing 2.5$ mm & $\varnothing 3.5$ mm Locking Screws
- Lateral, Midshaft & Hook Plates available

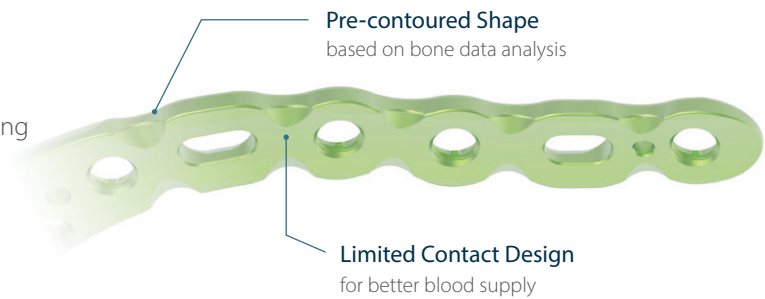


Material

Titanium 6Al 4V Alloy, ASTM F-136


Lateral Plate

- Fracture of the lateral clavicle

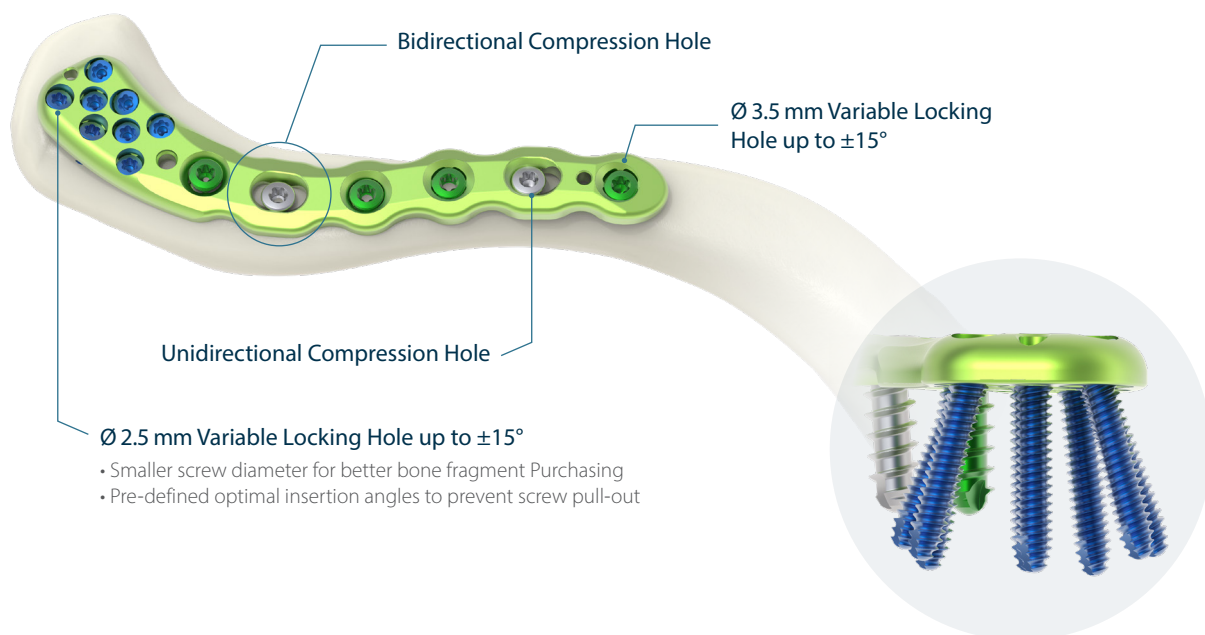

Midshaft Plate

- Fractures of the clavicle shaft
- Fracture of the clavicular midshaft

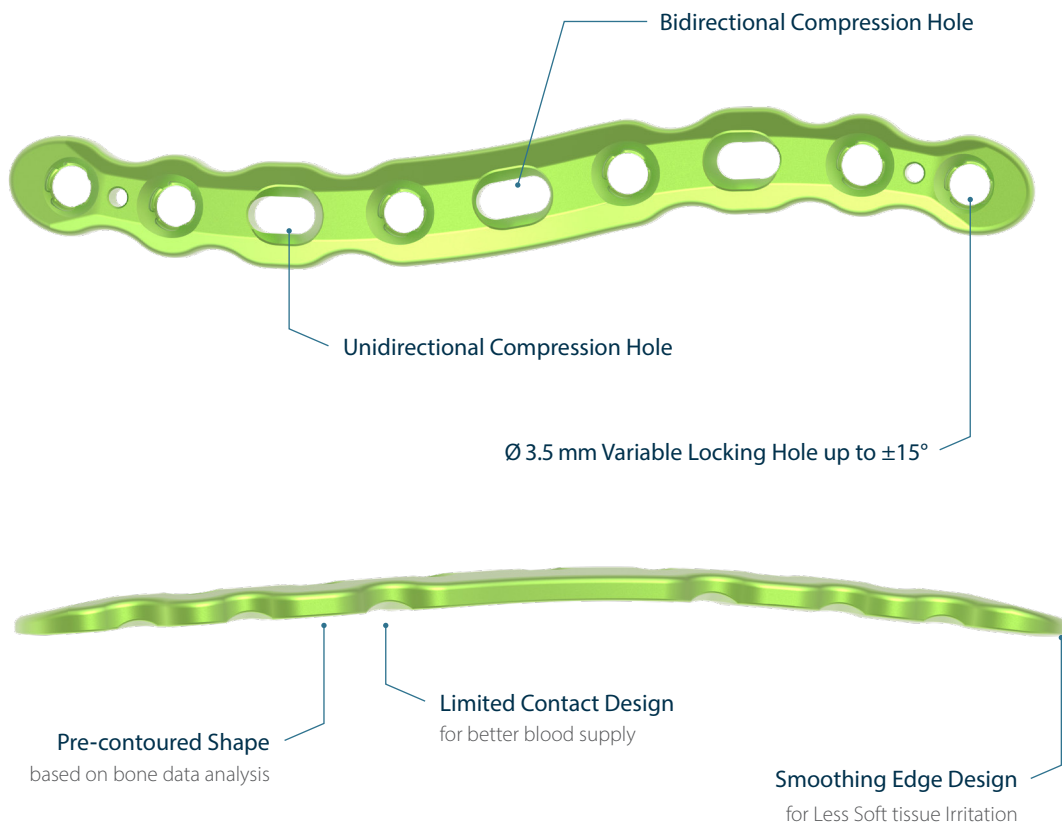



Two Variations of Curvatures for Better Fitting

- Increased Curvature
- Decreased Curvature

Features - continued


Hook Plate

- Acromioclavicular joint dislocation
- Fractures of the lateral clavicle

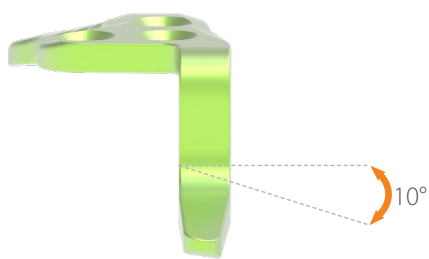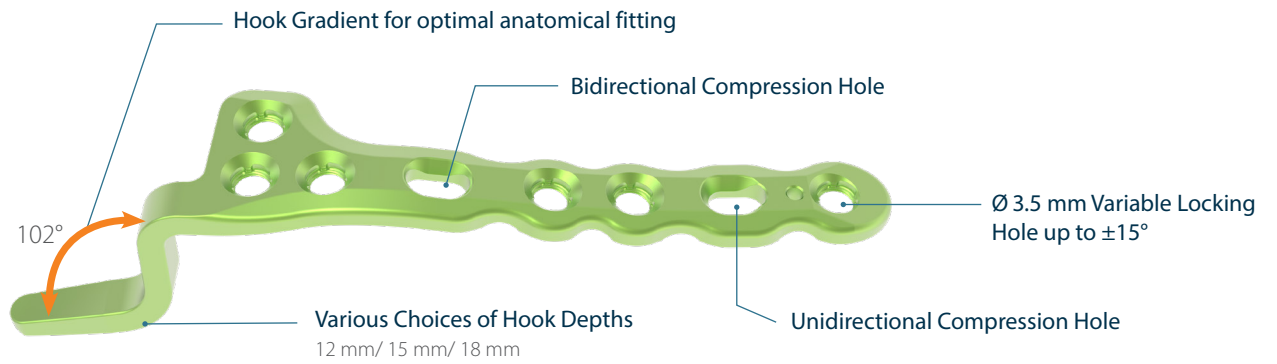


Hook Depths

- 12 mm
- 15 mm
- 18 mm

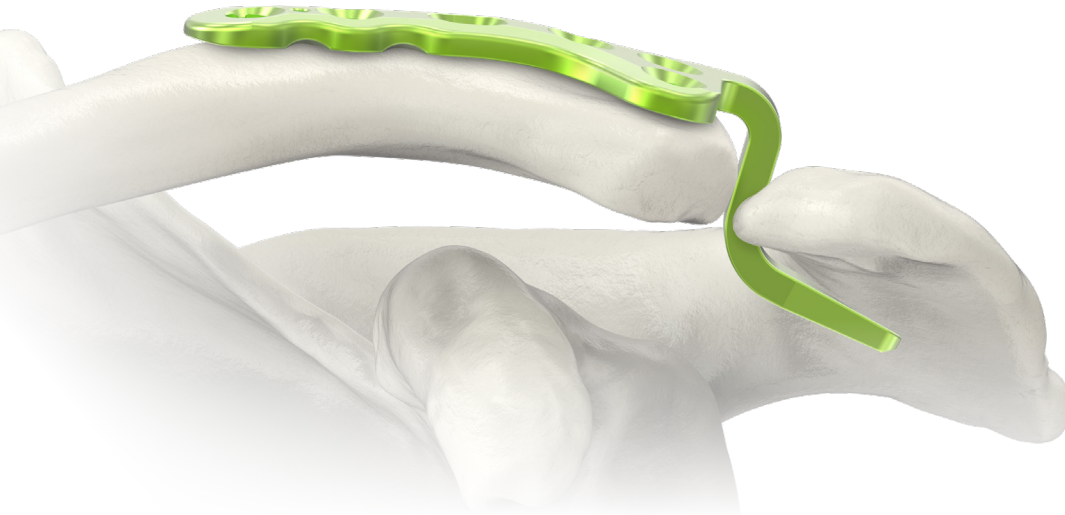

STARIX Screw with T10

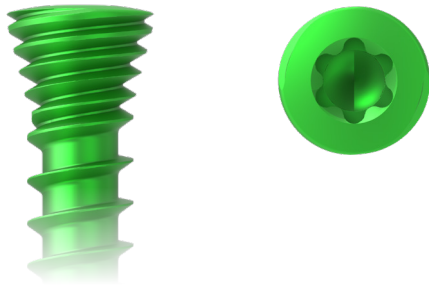

- Prevents Cam-out
- Allows Higher Torque Application
- Self-Retaining function
- Locking & Non-Locking Variations available


Variable Locking Interface with Locking Screw



- Poly-axial Screw Insertion
- Angle range: $\pm 15^\circ$
- Plate-Screw Locking Interface

Screw Options

with STARIX Pick-up

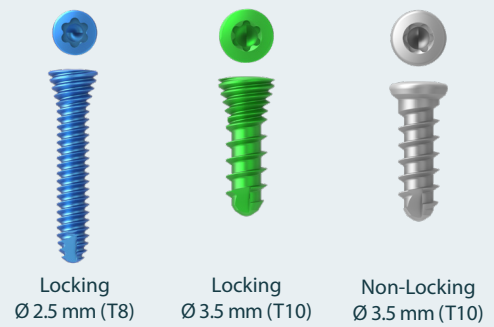


Locking
Ø 2.5 mm (T8)

Locking
Ø 3.5 mm (T10)

Non-Locking
Ø 3.5 mm (T10)


Drill Guide Block for Optimal Screw Insertion

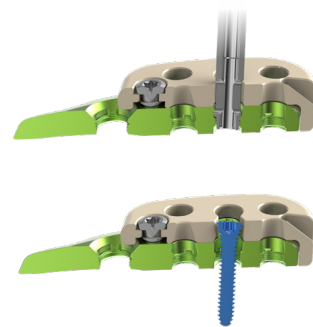

Drill Guide Block ensures optimal screw insertion angles


Clavicle Superior Lateral Plates

Left Plates

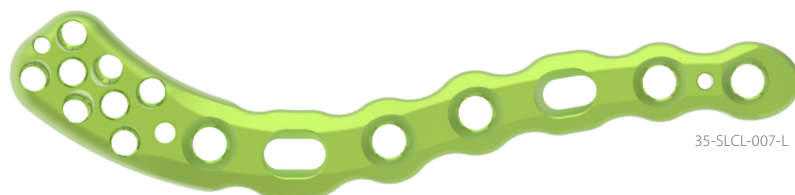

Code	Holes	Length	Thickness	Color	Set Quantity
35-SLCL-004-L	12	66.3 mm	3.5 mm	Light Green	2
35-SLCL-005-L	13	80.4 mm	3.5 mm	Light Green	2
35-SLCL-006-L	14	90.7 mm	3.5 mm	Light Green	2
35-SLCL-007-L	15	102.1 mm	3.5 mm	Light Green	2
35-SLCL-008-L	16	113.4 mm	3.5 mm	Light Green	2
35-SLCL-009-L	17	124.1 mm	3.5 mm	Light Green	2


Right Plates

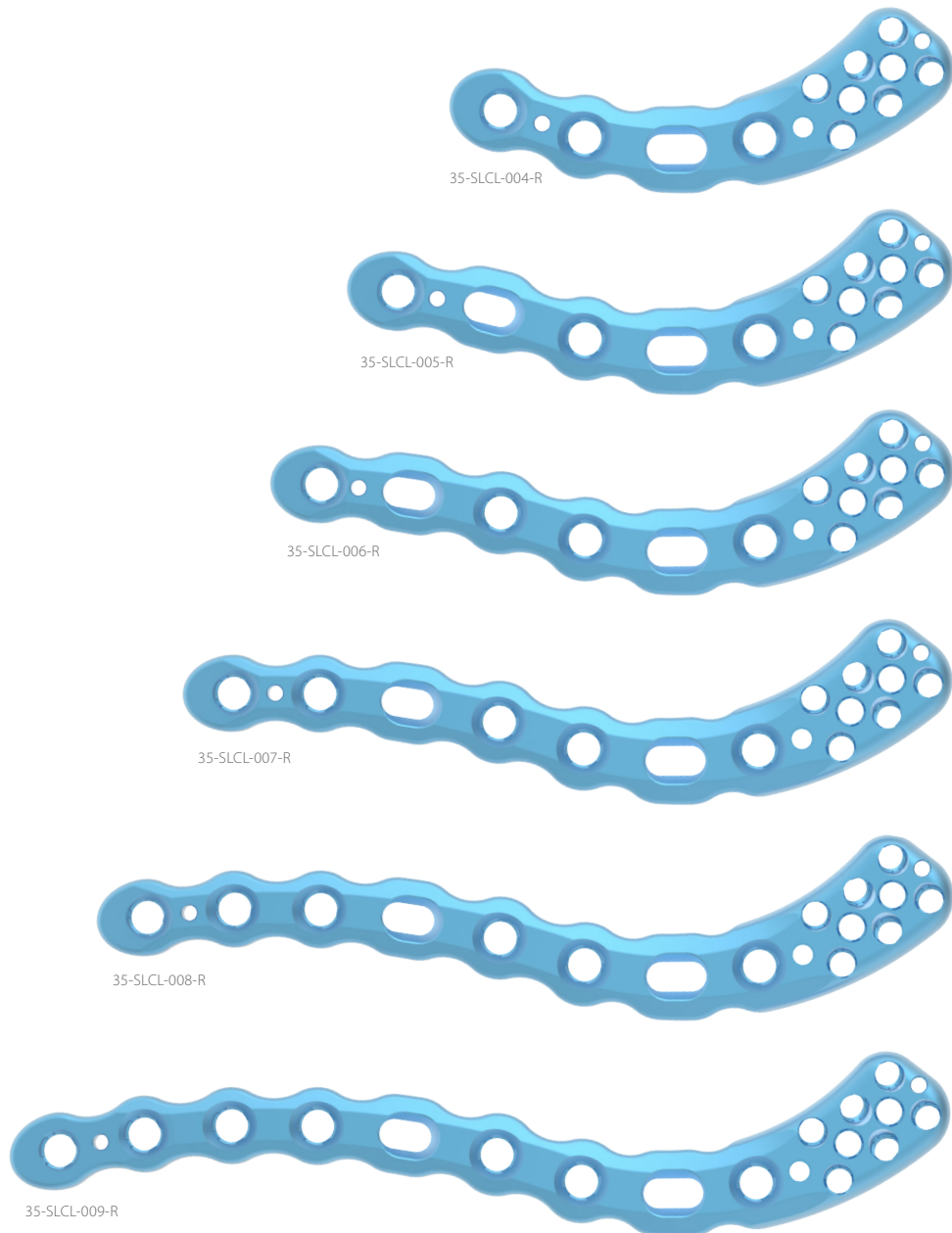

Code	Holes	Length	Thickness	Color	Set Quantity
35-SLCL-004-R	12	66.3 mm	3.5 mm	Light Blue	2
35-SLCL-005-R	13	80.4 mm	3.5 mm	Light Blue	2
35-SLCL-006-R	14	90.7 mm	3.5 mm	Light Blue	2
35-SLCL-007-R	15	102.1 mm	3.5 mm	Light Blue	2
35-SLCL-008-R	16	113.4 mm	3.5 mm	Light Blue	2
35-SLCL-009-R	17	124.1 mm	3.5 mm	Light Blue	2

Clavicle Superior Midshaft Plates

Decreased, Left Plates

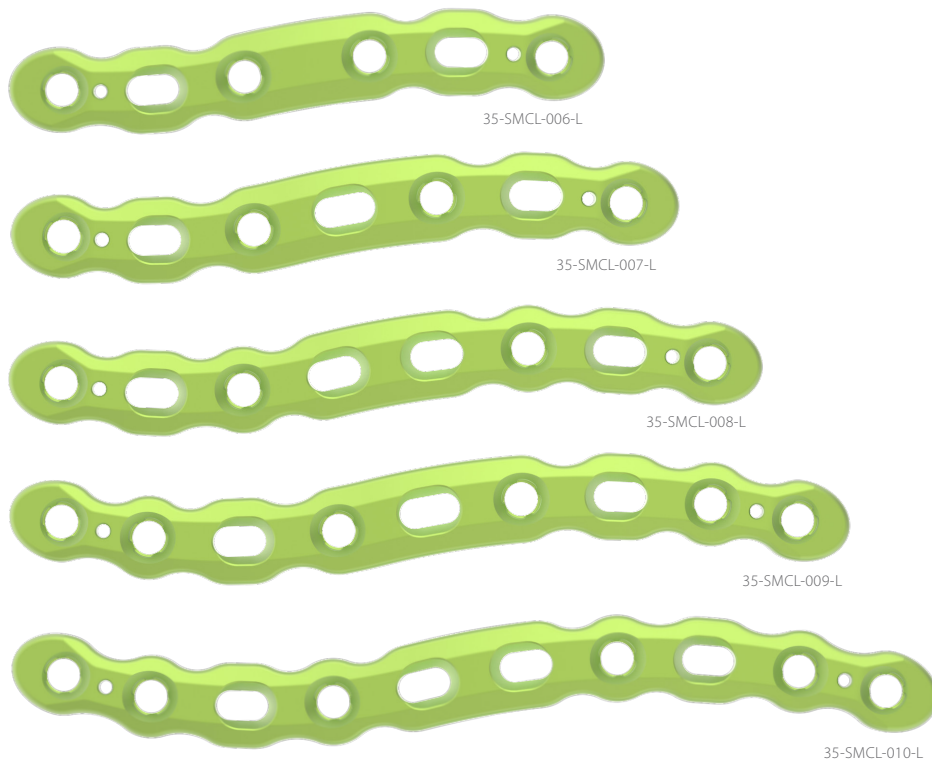


Code	Holes	Length	Thickness	Color	Set Quantity
35-SMCL-006-L	6	78.1 mm	3.5 mm	Light Green	2
35-SMCL-007-L	7	87.7 mm	3.5 mm	Light Green	2
35-SMCL-008-L	8	99.0 mm	3.5 mm	Light Green	2
35-SMCL-009-L	9	110.5 mm	3.5 mm	Light Green	2
35-SMCL-010-L	10	121.7 mm	3.5 mm	Light Green	2

Increased, Left Plates




Code	Holes	Length	Thickness	Color	Set Quantity
35-SMCL-008-LI	8	98.3 mm	3.5 mm	Light Green	2
35-SMCL-010-LI	10	120.5 mm	3.5 mm	Light Green	2

Decreased, Right Plates

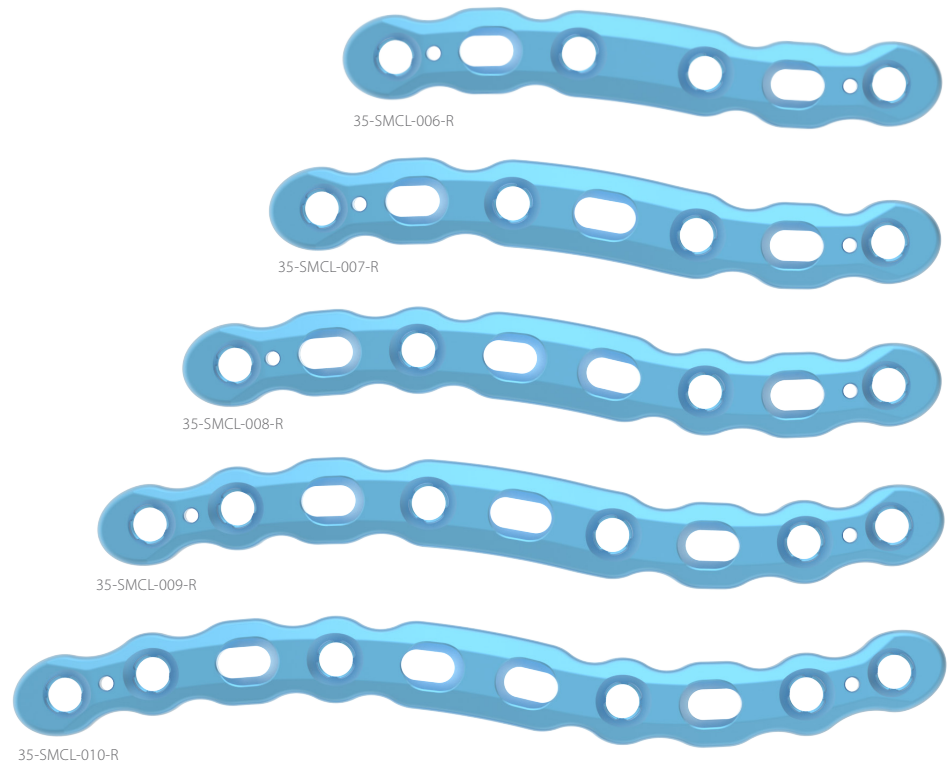

Code	Holes	Length	Thickness	Color	Set Quantity
35-SMCL-006-R	6	78.1 mm	3.5 mm	Light Blue	2
35-SMCL-007-R	7	87.7 mm	3.5 mm	Light Blue	2
35-SMCL-008-R	8	99.0 mm	3.5 mm	Light Blue	2
35-SMCL-009-R	9	110.5 mm	3.5 mm	Light Blue	2
35-SMCL-010-R	10	121.7 mm	3.5 mm	Light Blue	2


Increased, Right Plates

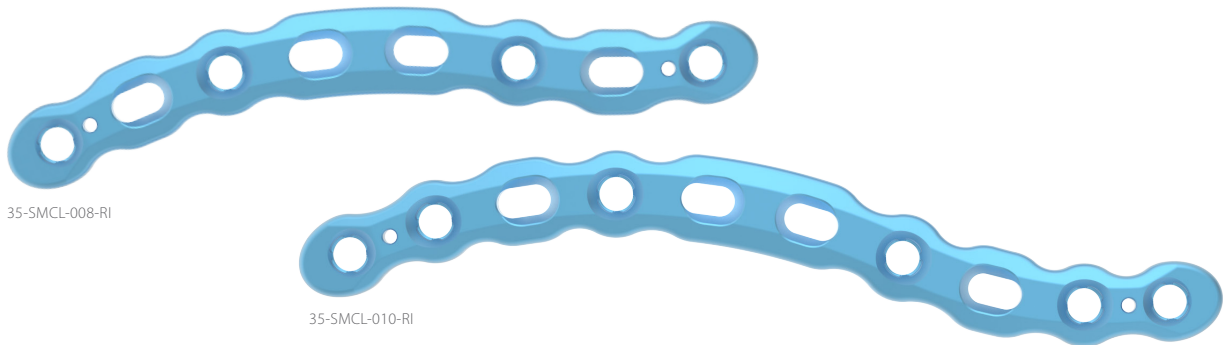

Code	Holes	Length	Thickness	Color	Set Quantity
35-SMCL-008-RI	8	98.3 mm	3.5 mm	Light Blue	2
35-SMCL-010-RI	10	120.5 mm	3.5 mm	Light Blue	2


Clavicle Hook Plates

Left Plates

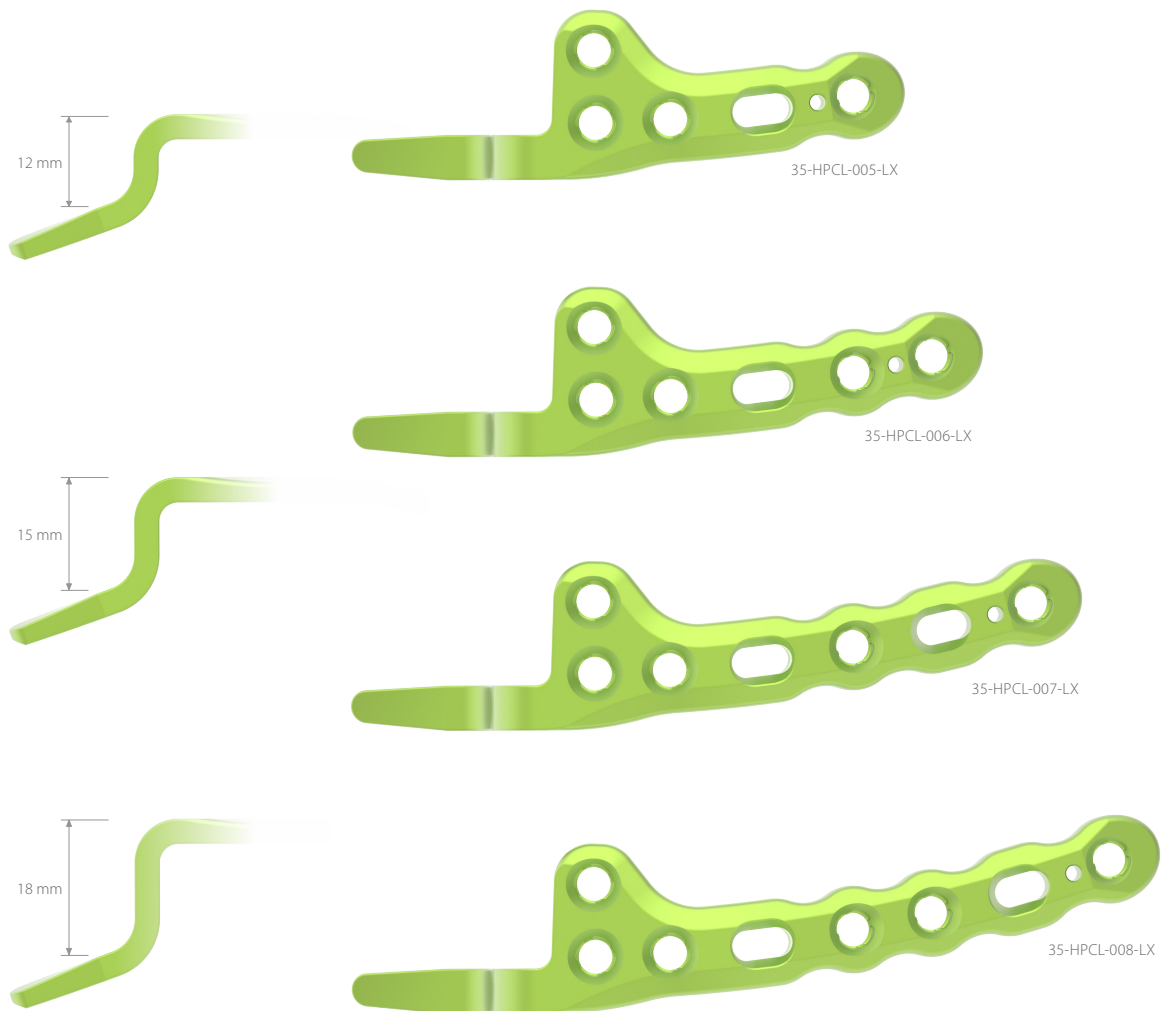

Code	Holes	Hook Depth	Length	Thickness	Color	Set Quantity
35-HPCL-005-L2	5	12 mm	72.7 mm	3.5 mm	Light Green	2
35-HPCL-006-L2	6	12 mm	82.9 mm	3.5 mm	Light Green	2
35-HPCL-007-L2	7	12 mm	96.0 mm	3.5 mm	Light Green	2
35-HPCL-008-L2	8	12 mm	105.1 mm	3.5 mm	Light Green	2
35-HPCL-005-L5	5	15 mm	72.7 mm	3.5 mm	Light Green	2
35-HPCL-006-L5	6	15 mm	82.9 mm	3.5 mm	Light Green	2
35-HPCL-007-L5	7	15 mm	96.0 mm	3.5 mm	Light Green	2
35-HPCL-008-L5	8	15 mm	105.1 mm	3.5 mm	Light Green	2
35-HPCL-005-L8	5	18 mm	72.7 mm	3.5 mm	Light Green	2
35-HPCL-006-L8	6	18 mm	82.9 mm	3.5 mm	Light Green	2
35-HPCL-007-L8	7	18 mm	96.0 mm	3.5 mm	Light Green	2
35-HPCL-008-L8	8	18 mm	105.1 mm	3.5 mm	Light Green	2

Right Plates

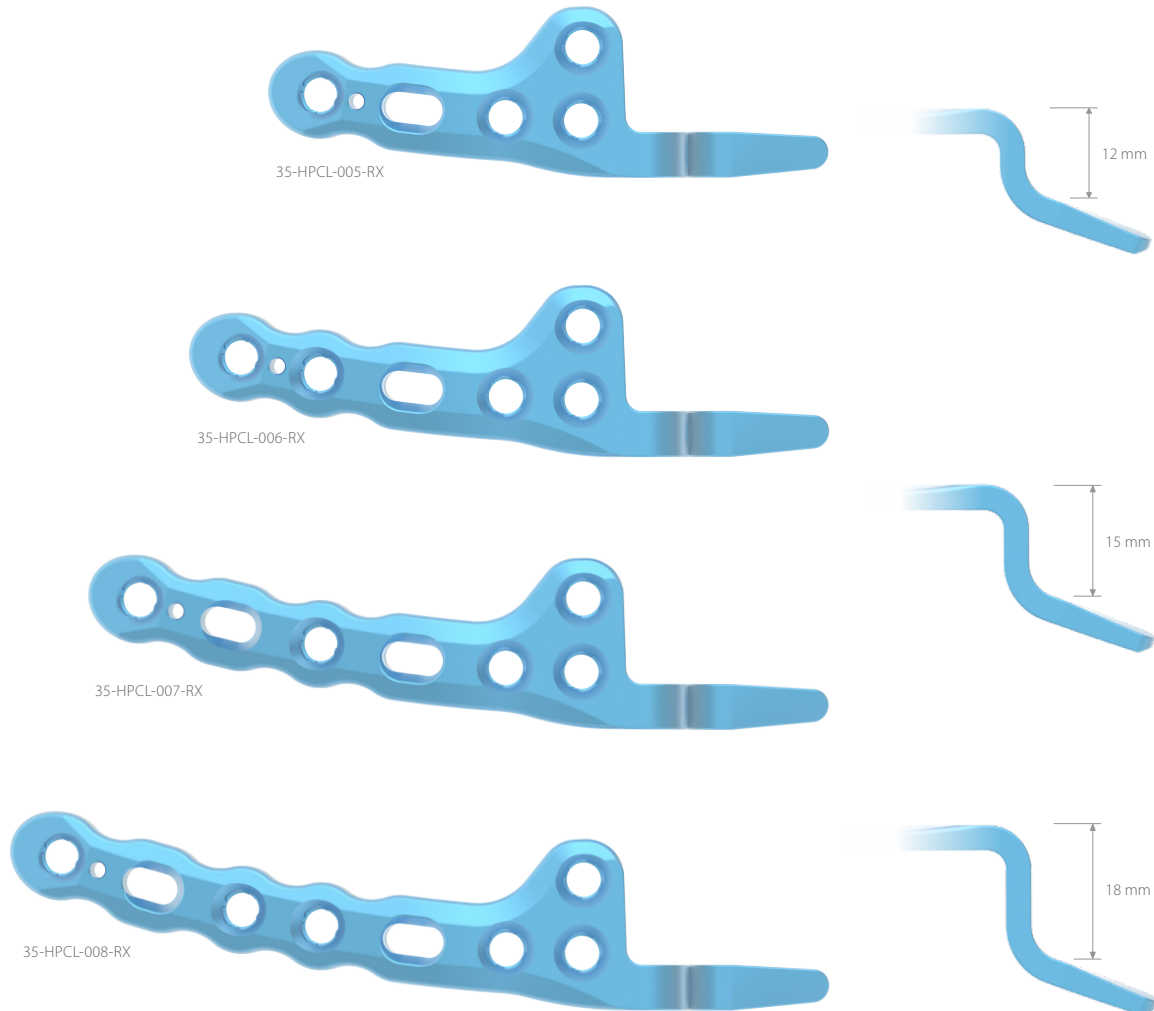

Code	Holes	Hook Depth	Length	Thickness	Color	Set Quantity
35-HPCL-005-R2	5	12 mm	72.7 mm	3.5 mm	Light Blue	2
35-HPCL-006-R2	6	12 mm	82.9 mm	3.5 mm	Light Blue	2
35-HPCL-007-R2	7	12 mm	96.0 mm	3.5 mm	Light Blue	2
35-HPCL-008-R2	8	12 mm	105.1 mm	3.5 mm	Light Blue	2
35-HPCL-005-R5	5	15 mm	72.7 mm	3.5 mm	Light Blue	2
35-HPCL-006-R5	6	15 mm	82.9 mm	3.5 mm	Light Blue	2
35-HPCL-007-R5	7	15 mm	96.0 mm	3.5 mm	Light Blue	2
35-HPCL-008-R5	8	15 mm	105.1 mm	3.5 mm	Light Blue	2
35-HPCL-005-R8	5	18 mm	72.7 mm	3.5 mm	Light Blue	2
35-HPCL-006-R8	6	18 mm	82.9 mm	3.5 mm	Light Blue	2
35-HPCL-007-R8	7	18 mm	96.0 mm	3.5 mm	Light Blue	2
35-HPCL-008-R8	8	18 mm	105.1 mm	3.5 mm	Light Blue	2

ARIX Clavicle System Screws

Ø 3.5 mm Locking Screws




Code	Length	Color	Set Quantity
35L-SO-L10-TA	10 mm	Green	10
35L-SO-L12-TA	12 mm	Green	10
35L-SO-L14-TA	14 mm	Green	10
35L-SO-L16-TA	16 mm	Green	10
35L-SO-L18-TA	18 mm	Green	10
35L-SO-L20-TA	20 mm	Green	10
35L-SO-L22-TA	22 mm	Green	10
35L-SO-L24-TA	24 mm	Green	10
35L-SO-L26-TA	26 mm	Green	10
35L-SO-L28-TA	28 mm	Green	10
35L-SO-L30-TA	30 mm	Green	10
35L-SO-L32-TA	32 mm	Green	10
35L-SO-L34-TA	34 mm	Green	10

Ø 3.5 mm Non-Locking Screws



Code	Length	Color	Set Quantity
35-SO-L10-T	10 mm	Silver	5
35-SO-L12-T	12 mm	Silver	5
35-SO-L14-T	14 mm	Silver	5
35-SO-L16-T	16 mm	Silver	5
35-SO-L18-T	18 mm	Silver	5
35-SO-L20-T	20 mm	Silver	5
35-SO-L22-T	22 mm	Silver	5
35-SO-L24-T	24 mm	Silver	5
35-SO-L26-T	26 mm	Silver	5
35-SO-L28-T	28 mm	Silver	5
35-SO-L30-T	30 mm	Silver	5
35-SO-L32-T	32 mm	Silver	5
35-SO-L34-T	34 mm	Silver	5

Ø 2.5 mm Locking Screws

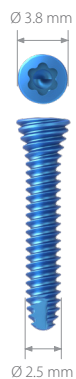

Code	Length	Color	Set Quantity
25L-SO-008-TA	8 mm	Blue	5
25L-SO-010-TA	10 mm	Blue	5
25L-SO-012-TA	12 mm	Blue	5
25L-SO-014-TA	14 mm	Blue	5
25L-SO-016-TA	16 mm	Blue	5
25L-SO-018-TA	18 mm	Blue	5
25L-SO-020-TA	20 mm	Blue	5
25L-SO-022-TA	22 mm	Blue	5
25L-SO-024-TA	24 mm	Blue	5
25L-SO-026-TA	26 mm	Blue	5

ARIX Clavicle System Instruments

Drill Sleeve & Drill Sleeve Handle



111-101



111-103



111-173



111-171



111-157

Code	Description	Set Quantity
111-101	Drill Sleeve for Ø 2.5 mm Screws, Fixed Angle	1
111-103	Drill Sleeve for Ø 2.5 mm Screws, Variable Angle	1
111-173	Drill Sleeve for Ø 3.5 mm Screws, Fixed Angle	1
111-171	Drill Sleeve for Ø 3.5 mm Screws, Variable Angle	1
111-157	Drill Sleeve Handle	2

Drill Bits



112-25-701



112-35-703

Code	Description	Set Quantity
112-25-701	Drill Bits for Ø 2.5 mm Screws(Drilling Length: 20 mm, Overall Length: 110 mm)	2
112-35-703	Drill Bits for Ø 3.5 mm Screws(Drilling Length: 25 mm, Overall Length: 118 mm)	2

Screwdriver Shaft



113-HF-613



113-HF-619

Code	Description	Set Quantity
113-HF-613	Screwdriver Shaft for STARIX T8	2
113-HF-619	Screwdriver Shaft for STARIX T10	2

ARIX Clavicle System Instruments

Depth Gauge



111-075



111-086

Code	Description	Set Quantity
111-075	Depth Gauge for Ø 2.5 mm Screws	1
111-086	Depth Gauge for Ø 3.5 mm Screws	1

Drill Guide



111-080



111-089

Code	Description	Set Quantity
111-080	Drill Guide for Ø 2.5 mm Screws, Fixed & Variable Angle	1
111-089	Drill Guide for Ø 3.5 mm Screws, Fixed & Variable Angle	1

Guide Pin Dispenser & Guide Pin



111-068-3



111-096

Code	Description	Set Quantity
111-068-3	Guide Pin Ø 1.6 mm	10
111-096	Guide Pin Dispenser	1

Screwdriver Handle



111-092

Code	Description	Set Quantity
111-092	Screwdriver Handle	2

Bone Clamp & Reduction Forceps



111-134



111-154

Code	Description	Set Quantity
111-134	Bone Clamp	2
111-154	Reduction Forceps, Large	2

Drill Guide Block




111-196-L



111-196-R

Code	Description	Set Quantity
111-196-L	Drill Guide Block, Left	1
111-196-R	Drill Guide Block, Right	1

Plate Bender

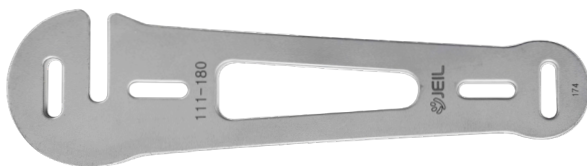

111-180

Code	Description	Set Quantity
111-180	Plate Bender	2

Clavicle Retractor




111-197

Code	Description	Set Quantity
111-197	Clavicle Retractor	1

ARIX Clavicle Set Configuration

Plates				Ø 3.5 mm Locking Screws			
No.	Code	Description	Set Quantity	No.	Code	Description	Set Quantity
1	35-SLCL-004-L	Clavicle Superior Lateral Plate, Left, 12 holes	2	1	35L-SO-L10-TA	Ø 3.5 mm Locking Screw, Length 10 mm	10
2	35-SLCL-005-L	Clavicle Superior Lateral Plate, Left, 13 holes	2	2	35L-SO-L12-TA	Ø 3.5 mm Locking Screw, Length 12 mm	10
3	35-SLCL-006-L	Clavicle Superior Lateral Plate, Left, 14 holes	2	3	35L-SO-L14-TA	Ø 3.5 mm Locking Screw, Length 14 mm	10
4	35-SLCL-007-L	Clavicle Superior Lateral Plate, Left, 15 holes	2	4	35L-SO-L16-TA	Ø 3.5 mm Locking Screw, Length 16 mm	10
5	35-SLCL-008-L	Clavicle Superior Lateral Plate, Left, 16 holes	2	5	35L-SO-L18-TA	Ø 3.5 mm Locking Screw, Length 18 mm	10
6	35-SLCL-009-L	Clavicle Superior Lateral Plate, Left, 17 holes	2	6	35L-SO-L20-TA	Ø 3.5 mm Locking Screw, Length 20 mm	10
7	35-SLCL-004-R	Clavicle Superior Lateral Plate, Right, 12 holes	2	7	35L-SO-L22-TA	Ø 3.5 mm Locking Screw, Length 22 mm	10
8	35-SLCL-005-R	Clavicle Superior Lateral Plate, Right, 13 holes	2	8	35L-SO-L24-TA	Ø 3.5 mm Locking Screw, Length 24 mm	10
9	35-SLCL-006-R	Clavicle Superior Lateral Plate, Right, 14 holes	2	9	35L-SO-L26-TA	Ø 3.5 mm Locking Screw, Length 26 mm	10
10	35-SLCL-007-R	Clavicle Superior Lateral Plate, Right, 15 holes	2	10	35L-SO-L28-TA	Ø 3.5 mm Locking Screw, Length 28 mm	10
11	35-SLCL-008-R	Clavicle Superior Lateral Plate, Right, 16 holes	2	11	35L-SO-L30-TA	Ø 3.5 mm Locking Screw, Length 30 mm	10
12	35-SLCL-009-R	Clavicle Superior Lateral Plate, Right, 17 holes	2	12	35L-SO-L32-TA	Ø 3.5 mm Locking Screw, Length 32 mm	10
13	35-SMCL-006-L	Clavicle Superior Midshaft Plate Decreased, Left, 6 holes	2	13	35L-SO-L34-TA	Ø 3.5 mm Locking Screw, Length 34 mm	10
14	35-SMCL-007-L	Clavicle Superior Midshaft Plate Decreased, Left, 7 holes	2				
15	35-SMCL-008-L	Clavicle Superior Midshaft Plate Decreased, Left, 8 holes	2	Ø 3.5 mm Non-Locking Screws			
16	35-SMCL-009-L	Clavicle Superior Midshaft Plate Decreased, Left, 9 holes	2	No.	Code	Description	Set Quantity
17	35-SMCL-010-L	Clavicle Superior Midshaft Plate Decreased, Left, 10 holes	2	1	35-SO-L10-T	Ø 3.5 mm Cortical Screw, Length 10 mm	5
18	35-SMCL-008-LI	Clavicle Superior Midshaft Plate Increased, Left, 8 holes	2	2	35-SO-L12-T	Ø 3.5 mm Cortical Screw, Length 12 mm	5
19	35-SMCL-010-LI	Clavicle Superior Midshaft Plate Increased, Left, 10 holes	2	3	35-SO-L14-T	Ø 3.5 mm Cortical Screw, Length 14 mm	5
20	35-SMCL-006-R	Clavicle Superior Midshaft Plate Decreased, Right, 6 holes	2	4	35-SO-L16-T	Ø 3.5 mm Cortical Screw, Length 16 mm	5
21	35-SMCL-007-R	Clavicle Superior Midshaft Plate Decreased, Right, 7 holes	2	5	35-SO-L18-T	Ø 3.5 mm Cortical Screw, Length 18 mm	5
22	35-SMCL-008-R	Clavicle Superior Midshaft Plate Decreased, Right, 8 holes	2	6	35-SO-L20-T	Ø 3.5 mm Cortical Screw, Length 20 mm	5
23	35-SMCL-009-R	Clavicle Superior Midshaft Plate Decreased, Right, 9 holes	2	7	35-SO-L22-T	Ø 3.5 mm Cortical Screw, Length 22 mm	5
24	35-SMCL-010-R	Clavicle Superior Midshaft Plate Decreased, Right, 10 holes	2	8	35-SO-L24-T	Ø 3.5 mm Cortical Screw, Length 24 mm	5
25	35-SMCL-008-RI	Clavicle Superior Midshaft Plate Increased, Right, 8 holes	2	9	35-SO-L26-T	Ø 3.5 mm Cortical Screw, Length 26 mm	5
26	35-SMCL-010-RI	Clavicle Superior Midshaft Plate Increased, Right, 10 holes	2	10	35-SO-L28-T	Ø 3.5 mm Cortical Screw, Length 28 mm	5
27	35-HPCL-005-L2	Clavicle Hook Plate, Left, Depth 12 mm, 5 holes	2	11	35-SO-L30-T	Ø 3.5 mm Cortical Screw, Length 30 mm	5
28	35-HPCL-006-L2	Clavicle Hook Plate, Left, Depth 12 mm, 6 holes	2	12	35-SO-L32-T	Ø 3.5 mm Cortical Screw, Length 32 mm	5
29	35-HPCL-007-L2	Clavicle Hook Plate, Left, Depth 12 mm, 7 holes	2	13	35-SO-L34-T	Ø 3.5 mm Cortical Screw, Length 34 mm	5
30	35-HPCL-008-L2	Clavicle Hook Plate, Left, Depth 12 mm, 8 holes	2				
31	35-HPCL-005-L5	Clavicle Hook Plate, Left, Depth 15 mm, 5 holes	2	Ø 2.5 mm Locking Screws			
32	35-HPCL-006-L5	Clavicle Hook Plate, Left, Depth 15 mm, 6 holes	2	No.	Code	Description	Set Quantity
33	35-HPCL-007-L5	Clavicle Hook Plate, Left, Depth 15 mm, 7 holes	2	1	25L-SO-008-TA	Ø 2.5 mm Locking Screw, Length 8 mm	5
34	35-HPCL-008-L5	Clavicle Hook Plate, Left, Depth 15 mm, 8 holes	2	2	25L-SO-010-TA	Ø 2.5 mm Locking Screw, Length 10 mm	5
35	35-HPCL-005-L8	Clavicle Hook Plate, Left, Depth 18 mm, 5 holes	2	3	25L-SO-012-TA	Ø 2.5 mm Locking Screw, Length 12 mm	5
36	35-HPCL-006-L8	Clavicle Hook Plate, Left, Depth 18 mm, 6 holes	2	4	25L-SO-014-TA	Ø 2.5 mm Locking Screw, Length 14 mm	5
37	35-HPCL-007-L8	Clavicle Hook Plate, Left, Depth 18 mm, 7 holes	2	5	25L-SO-016-TA	Ø 2.5 mm Locking Screw, Length 16 mm	5
38	35-HPCL-008-L8	Clavicle Hook Plate, Left, Depth 18 mm, 8 holes	2	6	25L-SO-018-TA	Ø 2.5 mm Locking Screw, Length 18 mm	5
39	35-HPCL-005-R2	Clavicle Hook Plate, Right, Depth 12mm, 5 holes	2	7	25L-SO-020-TA	Ø 2.5 mm Locking Screw, Length 20 mm	5
40	35-HPCL-006-R2	Clavicle Hook Plate, Right, Depth 12mm, 6 holes	2	8	25L-SO-022-TA	Ø 2.5 mm Locking Screw, Length 22 mm	5
41	35-HPCL-007-R2	Clavicle Hook Plate, Right, Depth 12mm, 7 holes	2	9	25L-SO-024-TA	Ø 2.5 mm Locking Screw, Length 24 mm	5
42	35-HPCL-008-R2	Clavicle Hook Plate, Right, Depth 12mm, 8 holes	2	10	25L-SO-026-TA	Ø 2.5 mm Locking Screw, Length 26 mm	5
43	35-HPCL-005-R5	Clavicle Hook Plate, Right, Depth 15mm, 5 holes	2				
44	35-HPCL-006-R5	Clavicle Hook Plate, Right, Depth 15mm, 6 holes	2				
45	35-HPCL-007-R5	Clavicle Hook Plate, Right, Depth 15mm, 7 holes	2				
46	35-HPCL-008-R5	Clavicle Hook Plate, Right, Depth 15mm, 8 holes	2				
47	35-HPCL-005-R8	Clavicle Hook Plate, Right, Depth 18mm, 5 holes	2				
48	35-HPCL-006-R8	Clavicle Hook Plate, Right, Depth 18mm, 6 holes	2				
49	35-HPCL-007-R8	Clavicle Hook Plate, Right, Depth 18mm, 7 holes	2				
50	35-HPCL-008-R8	Clavicle Hook Plate, Right, Depth 18mm, 8 holes	2				

Instruments			
No.	Code	Description	Set Quantity
1	112-108	ARIX Clavicle System Tray	1
2	111-101	Drill Sleeve for Ø 2.5 mm Screws, Fixed Angle	1
3	111-103	Drill Sleeve for Ø 2.5 mm Screws, Variable Angle	1
4	111-173	Drill Sleeve for Ø 3.5 mm Screws, Fixed Angle	1
5	111-171	Drill Sleeve for Ø 3.5 mm Screws, Variable Angle	1
6	111-157	Drill Sleeve Handle	2
7	112-25-701	Drill Bits for Ø 2.5 mm Screws	2
8	112-35-703	Drill Bits for Ø 3.5 mm Screws	2
9	113-HF-613	Screwdriver Shaft for STARIX T8	2
10	113-HF-619	Screwdriver Shaft for STARIX T10	2
11	111-075	Depth Gauge for Ø 2.5 mm Screws	1
12	111-086	Depth Gauge for Ø 3.5 mm Screws	1
13	111-080	Drill Guide for Ø 2.5 mm Screws, Fixed & Variable Angle	1
14	111-089	Drill Guide for Ø 3.5 mm Screws, Fixed & Variable Angle	1
15	111-068-3	Guide Pin Ø 1.6 mm	10
16	111-096	Guide Pin Dispenser	1
17	111-092	Screwdriver Handle	2
18	111-134	Bone Clamp	2
19	111-154	Reduction Forceps, Large	2
20	111-196-L	Drill Guide Block Left	1
21	111-196-R	Drill Guide Block Right	1
22	111-180	Plate Bender	2
23	111-197	Clavicle Retractor	1

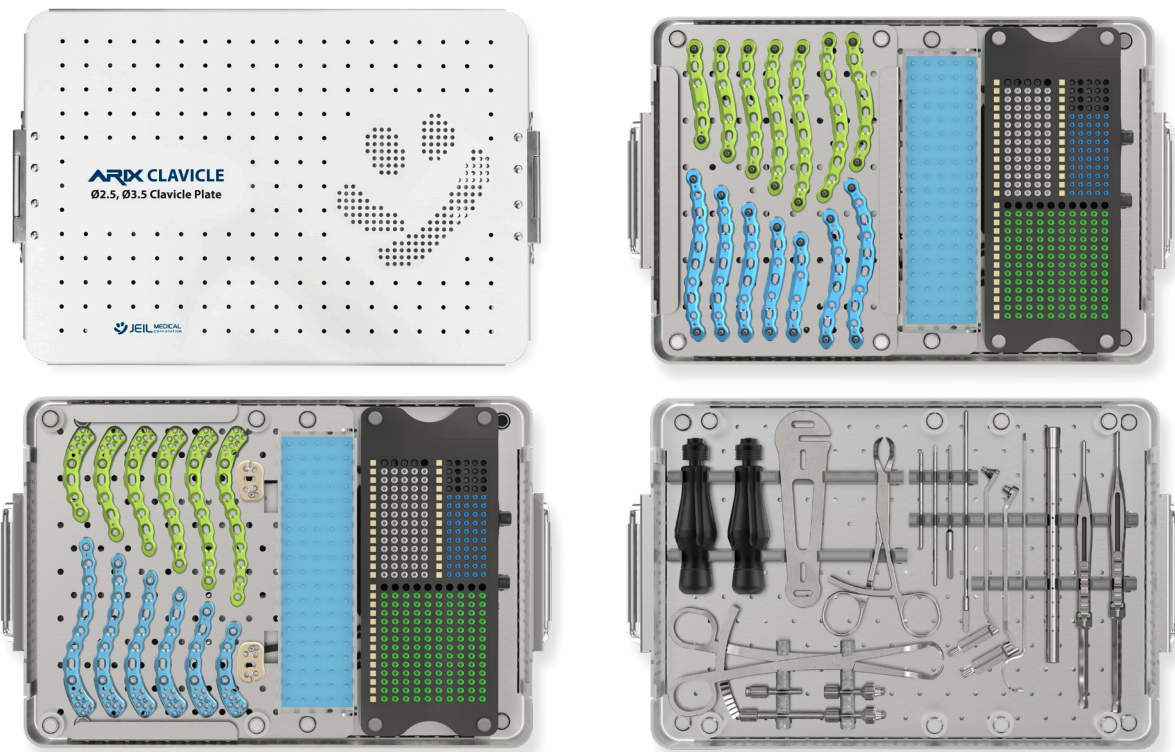

112-108
ARIX Clavicle System Tray

JEIL MEDICAL CORPORATION

702•703•704•705•706•804•805•807•812•815-ho, 55, Digital-ro 34-gil, Guro-gu, Seoul, 08378, Rep. of Korea

DOMESTIC SALES_ **T** 82.2.850.3500 **F** 82.2.850.3535 **E** jeil@jeilmed.co.kr
INTERNATIONAL SALES_ **T** 82.2.850.3935 **F** 82.2.850.3537 **E** global@jeilmed.co.kr

www.jeilmed.co.kr